




## **Fixed Height Aluminium Bed Cradle**

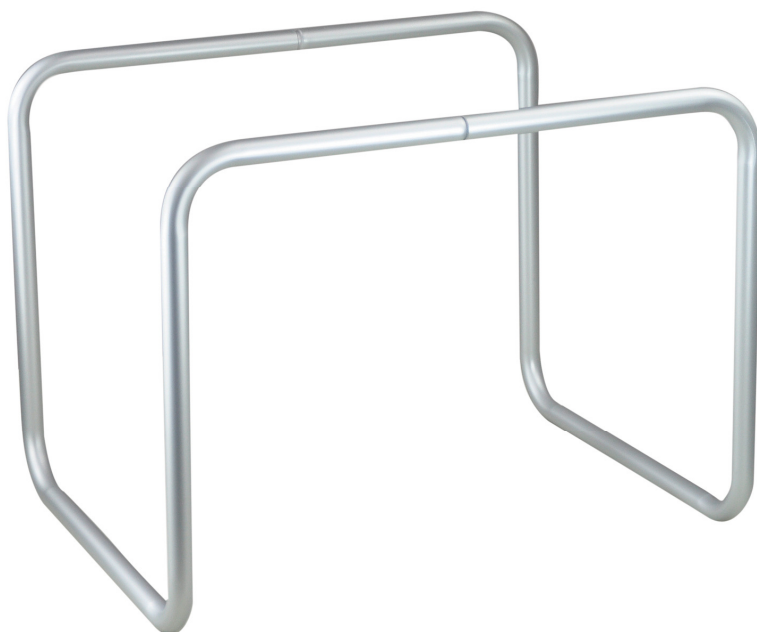

**VY441**

## **Assembly and Operating Instructions**

**The Fixed Height Aluminium Bed Cradle from Aidapt is a metal frame designed to keep bed covers away from the body. It is useful for those suffering from ulcers or skin irritation.**

## **INTRODUCTION**

Thank you for deciding to purchase the Fixed Height Aluminium Bed Cradle from Aidapt. The Bed Cradle is manufactured from the finest material and components available, when installed and used correctly it is designed to give many years of reliable trouble-free service. The Bed Cradle can be fitted to single, double, queen and king size beds.

## **INSTALLATION**

The base is worked between the mattress and bed base and can be used at the side or end of the bed.

## **IMPORTANT**

The Bed Cradle is not intended as an aid to help the user in and out of bed.

***NB. Do not lean on the Bed Cradle or apply heavy loads.***

The Bed Cradle has a **maximum load of 5 kg** and is designed to hold bed clothes only.


## **CLEANING**

The Adjustable Bed Cradle should be cleaned using a non-abrasive cleaner or mild detergent with a soft cloth only. Abrasive cleaners e.g. AJAX and/ or abrasive cleaning pads will seriously damage the product beyond repair and should not be used. Always ensure that you wipe dry the equipment after cleaning.

## **MAINTENANCE**

Please check the cradle periodically for any signs of damage and if in doubt contact your supplier.

## TECHNICAL SPECIFICATION

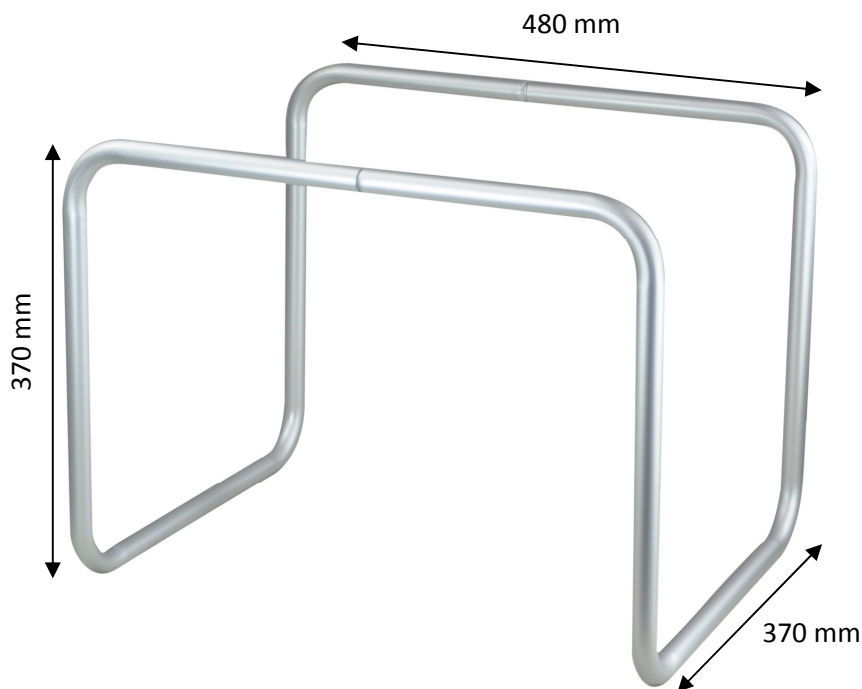

## IMPORTANT INFORMATION

The information given in this instruction booklet must not be taken as forming part of or establishing any contractual or other commitment by Aidapt Bathrooms Limited or Aidapt (Wales) Ltd or its agents or its subsidiaries and no warranty or representation concerning the information is given.

Please exercise common sense and do not take any unnecessary risks when using this product; as the user you must accept liability for safety when using the product.

**Aidapt Bathrooms Ltd, Lancots Lane, Sutton Oak, St Helens, WA9 3EX**  
**Telephone: +44 (0) 1744 745 020 • Fax: +44 (0) 1744 745 001 • Web: [www.aidapt.co.uk](http://www.aidapt.co.uk)**  
**Email: [accounts@aidapt.co.uk](mailto:accounts@aidapt.co.uk) • [adaptations@aidapt.co.uk](mailto:adaptations@aidapt.co.uk) • [sales@aidapt.co.uk](mailto:sales@aidapt.co.uk)**